Cedis[®]

Ultra Clean

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
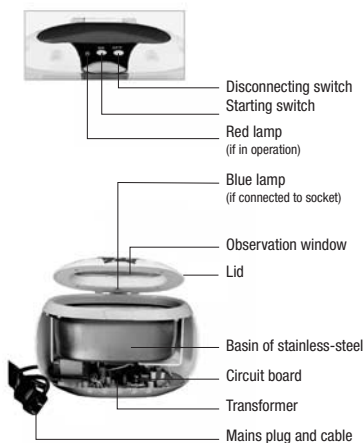
Before first use:

- Note mains voltage of device.
- Place device for operation on a dry and solid surface.
- Do not dip device in liquid.
- During operation of ultrasonic device do not put your hand into basin.
- Only meant for cleaning earmolds, glasses and jewelry – not meant for cleaning hearing aids or hearing aid glasses!
- Do not use inflammable liquids in cleansing basin.
- Do not pull out mains plug with cable out of the socket.
- In case of damaged cable immediate exchange is necessary by specialist.
- Before emptying the device, it is absolutely necessary to pull the mains plug.
- Keep away from heat sources.
- Do not force any objects through openings.
- Keep away from children.

Instruction for the Cleaning of Earmolds

1. Open the lid of the device and take out the insert. Fill the device with luke-warm water until the marking and add one measuring spoon of the special cleansing powder for the cleansing of earmolds.

ATTENTION:
Never put device
into operation without water!

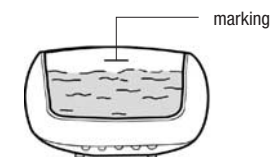
2. Separate the earmold from the hearing aid. Put the earmold into the insert and together into the device. ATTENTION: Never put the hearing aid into the cleansing liquid!
 3. Close the lid. Plug mains into the socket. The blue basin lamp is on.
 4. Press the ON-switch to put the device into operation (the red lamp is on).
 5. After 3 minutes the device switches off automatically. To repeat the cleansing process (in case of heavy soiling) press the ON-switch once again. At any time, you can end the cleansing process yourself by pressing the OFF-switch.
 6. After the cleansing process has come to an end, take out the earmold and rinse it with running water.
 7. Dry the earmold and blow through the venting with an air blower if necessary. There should not be any residual liquid in the tube or venting.
 8. Join the earmold to the hearing aid.
- ATTENTION: Do not use the cleansing water again!

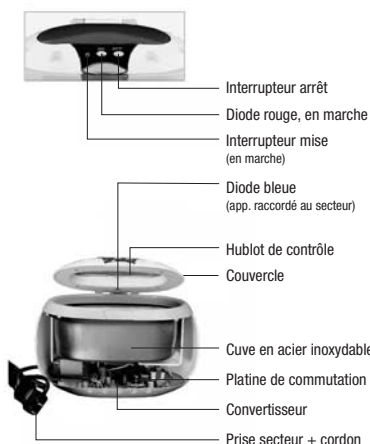For the cleaning of glasses and jewelry please use cleansing liquids from a specialist shop or ask your optician or at the jeweler's.

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egger Otoplastik + Labortechnik GmbH · Aybühlweg 59 · D-87439 Kempten/Germany